Set n' Go Height Adjustable Rollators Item RJ4700, RJ4718 (WIDE)

The Lumex Set n' Go Height Adjustable Rollator offers the convenience and flexibility of a two-in-one product all rolled into one. This accommodates for a broad range of user heights which in turn offers more comfort for the user when seated.

- Easy height adjustment on the frame allows seat height to be adjusted from 18" to 22"(19.5" to 23.5" for WIDE model) in 1" increments
- Comfortable, wide, seamless padded seat with zippered pouch under seat for security of personal items (15^{3/4"} on WIDE)
- Also includes black pouch for storage of larger items
- Adjustable handle height with knobs to secure in place
- Easy to operate locking and "quick stop" loop brakes
- 6" wheel (8" on WIDE) • 300 lb maximum weight capacity, (350 lb on WIDE),
- EVENLY DISTRIBUTED • Designed to meet the requirements of HCPCS: E0143 and E0156

Sure-GaitTMII Three Wheeled Steel Rollator Item 609101A

As a steel frame, three wheel rollator, the Sure-Gait II provides maximum maneuverability and comfort for those who do not require the extra stability of a four wheel rollator.

- Chrome-plated steel frame
- Adjustable handle height
- Easy one handed folding design
- Locking loop brakes
- Solid, puncture resistant 8" wheels for indoor and outdoor use
- Carry pouch included
- 300 lb maximum weight capacity EVENLY DISTRIBUTED
- Designed to meet the requirements of HCPCS code E0143

Information contained herein is subject to change. The most current and complete product information can be found on our website. www.grahamfield.com

© 2014, GF Health Products, Inc. All Rights Reserved. Lumex is a trademark of GF Health Products, Inc. GF Health Products, Inc. is an ISO 13485 Certified Company.

GF 3-Wheel Cruiser Item 609201

Bondi Blue

Metallic Burgundy

Peweter Metallic

The GF 3-Wheel Cruiser is designed to provide convenient mobility for those individuals who do not require a seat or the additional stability of a four-wheel rollator. This new light weight aluminum walking aid allows the user to have more freedom in movement to perform daily activities.

- Lightweight aluminum frame weighs only 10.5 lbs
- Ergonomic hand grips with locking loop brakes
- Adjustable handle height for broad range of users
- 8" black wheels for indoor and outdoor use
- Comes complete with zippered carry pouch
- 250 lb maximum weight capacity EVENLY DISTRIBUTED
- Available in 3 new attractive colors: Pewter Metallic, Metallic Burgundy and Bondi Blue
- Designed to meet the requirements of HCPCS code E0143

Hybrid Rollator Transport Chair Item LX1000

An innovative rollator that combines the function of a wheeled walker and a transport chair into one lightweight design

- Lightweight frame durable, aluminum frame weighs only 22 lbs (this includes footrests) making it easier for users to move about freely
- Hybrid concept designed to be used primarily as a rollator while easily transforming into a transport chair with removable footrests (adjustable in length) when user needs a short rest
- Ultra-Compact design tool-free assembly and folds effortlessly for convenient transport and storage
- Easy-to-use handbrakes allow for short and longer period braking when used as rollator and transport chair
- Flexible, soft polyester backrest and 18" wide mesh seat with pull up strap for easy folding provides comfort for user
- Large 8" wheels for easy rolling on most surfaces • Complete with storage
- pouch for personal items Available in attractive, contemporary colors: Titanium and Majestic Blue • Designed to meet the requirements of HCPCS E0143 & E0156



Get What You Need... & A Whole Lot More From Lumex[®], The Most Trusted Name In Walking Aids.



GF1000097 RevE16

graham field



The Lumex WalkaboutTM Four-Wheel Rollators provide the ideal combination of comfort, appeal, and affordability. The Walkabout series offers both high quality and a variety of options to meet the needs of a wide spectrum of patients. No two patients are alike, so let Lumex, with our complete line of rollators, help put their lives back in motion.

Lumex[®] Rollator Features

tem RJ4400

Designed with an aluminum frame built for durability, these lightweight rollators allow for easy maneuverability and fold quickly and effortlessly for transport. Each rollator features a wide padded seat and tool-free removable padded backrest for enhanced user comfort. All items have ergonomic hand grips that provide a better grip as well as locking loop brakes that are easy to operate and locking rear wheels for stable positioning when sitting.

Item RJ4405

The adjustable handle height accommodates a broad range of user heights with hemi models available for shorter individuals. Each comes with a wire basket that can be placed underneath the seat or, on some models, in front for carrying personal items. The Walkabout series is designed to meet the requirements of HCPCS codes E0143 and E0156 (E0149 for Bariatric items) and offers a Limited Lifetime Warranty as stated on the laminated color operating instructions included.

Walkabout[™] Imperial Bariatric Items RJ4400, RJ4405

The Walkabout Imperial rollator is designed for the bariatric user, featuring a 20"-wide padded seat and offers a 450 lb maximum weight capacity EVENLY DISTRIBUTED. It also features heavy-duty 8" wheels that are ideal for both indoor and outdoor use. Offered in both standard and ConTour options.

Walkabout[™] Deluxe ConTour Item RJ4805

The Walkabout Deluxe ConTour rollators offer the same benefits as the Walkabout Lite rollator, while featuring a contoured padded backrest and heavy-duty 8" wheels that are ideal for both indoor and outdoor use. Also available is an add on adjustable backrest (RJ4805-ABAR) that allows the user to adjust the backrest for ideal positioning.

Walkabout[™]Lite Item RJ4300

Our most popular item weighs only 14.5 lbs and is offered in various attractive colors: teal green, burgundy, pink, royal blue, aqua, lavender, champagne, and black. Each has 6" wheels and has a 300 lb. maximum weight capacity EVENLY DISTRIBUTED.

Item RJ4300


The walkabout lite folds effortlessly with folding strap for convenient transport and storage.

Walkabout Wide Item RJ4318

Designed for the user who is "in between sizes" of a standard or bariatric rollator, this Walkabout Wide offers more room and added comfort. The wider frame features an 18.5" wide opening in between the handgrips, a 15.5" wide padded seat and a 350 lb maximum weight capacity EVENLY DISTRIBUTED.



WIDE rollators offer 18.5" between handgrips and larger seats for added comfort

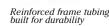

Item RJ4805



Item RJ4302

Easy grip knob makes for simple backbar removal for storage and transport.

Available colors *Color offerings vary for different models





As part of our Quality initiative, our Walkabout series meets ISO standards.



Walkabout[™]Lite Hemi Items RJ4302, RJ44002

This Hemi rollator offers the same benefits as the Walkabout Lite rollator yet is designed with a lower seat to floor height to accommodate shorter individuals. 6" wheels are available on this item. An Imperial Bariatric Hemi is also available (RJ4402R) to support a 450 lb. maximum weight capacity EVENLY DISTRIBUTED.

Walkabout[™]Junior Item RJ4301

This lightweight Walkabout Junior rollator is ideal for the petite or junior user. It is designed with a smaller frame, weighs only 12.5 lbs and has 6" wheels.

Aqua

Lavender

Champagne

Black



Rear Reflectors increase visibility at night!