

INVACARE®
NUTRON® R51™
POWER WHEELCHAIR



For a highly maneuverable power chair in an attractive package, look no further than the **Invacare® Nutron® R51™ and R51™LXP Power Wheelchairs**. Each of these models feature a rear-wheel drive wheelbase, enabling them to easily navigate obstacles inside and outside the home. Each chair folds for easy transport when the consumer's active lifestyle takes him or her on the road. The LXP package increases the motor output to 650 watts for a higher maximum speed and a weight capacity of 300 lb.

Features

Traditional rear-wheel drive and compact size provide excellent maneuverability and accessibility

Oversize, high-strength tubular frame with bold styling

Seat-to-floor heights as low as 17.75" permit easier transfers

Reinforced crossbrace design folds to 13° for transport

Back angle and height adjust for maximum comfort and support


Telescoping front riggings allow optimal foot placement

Flat-free wheel inserts mean worry-free operation

Powerful 350 or 650-watt, high-torque motors ensure dependable performance

Regenerative braking and automatic motor locks for smooth stopping

HCPCS Code: K0822

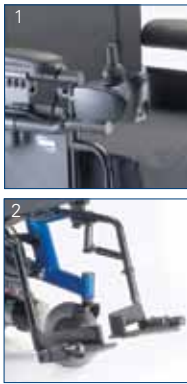

1. Invacare® MK5™ Electronics
2. Telescoping front riggings



Invacare Corporation
www.invacare.com

USA

One Invacare Way
Elyria, Ohio
44035-4190
(800) 333-6900

Canada

570 Matheson Blvd. E.,
Unit 8
Mississauga, Ontario
L4Z 4G4 Canada
(800) 668-5324

©2012

Invacare Corporation.
All rights reserved.
Trademarks are identified
by the symbols ™ and
®. All trademarks are
owned by or licensed
to Invacare Corporation
unless otherwise noted.
Specifications are subject
to change without
notification.
Form No. 01-337
Rev. 10/12

Model #	Overall Dimensions	Seat Dimensions	Seat-to-Floor Height	Drive Wheels	Ground Clearance	Turning Radius	Incline
R51, R51LXP	41.3" L 25" W 35.2" H	14"-20" W 16"-18" D	17.75"-19.75"	12.5" x 2.25"	3.5"	32.9"	9°
	Total Weight	Weight Capacity	Motor	Maximum Speed	Battery	Battery Range	Limited Warranty
	169 lb. 186 lb. (LXP)	250 lb. 300 lb. (LXP)	2 Pole 4 Pole (LXP)	4 mph 6 mph (LXP)	U1 22NF	N/A	60 Months Base Frame 36 Months Seat Frame 12 Months Electronics/ Motor/Gearbox

All references to HCPCS codes: Providers are responsible for determining the appropriate billing codes when submitting claims to the Medicare program, and should consult an attorney or other advisor to discuss specific situations in further detail.

